

Memo-Form Plus

AHSM-069 Rev 2.0 (24.09.24)

Static Replacement Mattress



The **Memo-Form Plus** mattress has been designed for users considered to be at 'High Risk' of developing pressure ulcers, catering for both acute healthcare and community environments.

Comprising of a visco elastic castellated foam insert, which allows both immersion and support, the **Memo-Foam Plus** mattress effectively re-distributes pressure whilst offering high levels of comfort and stability. The temperature sensitive foam conforms to an individual's shape, creating a greater body contact area, which leads to effective pressure redistribution.


- Castellated visco foam contributes to an effective pressure ulcer prevention plan, offering immersion, comfort, and postural support
- Internal air channels in the foam's surface help to reduce the buildup of heat and moisture on the skin
- Smooth multi-stretch PU cover is both waterproof and vapour permeable, reducing shear and friction
- High frequency welded seams reduce the risk of fluid ingress, supporting infection control protocols
- The non-turn mattress system reduces manual handling risks
- Reinforced side wall and base offer stability and support which assists with patient transfers
- Suitable for electric profiling bed frames
- Available in single and double sizes, bespoke sizes are also available upon request

Castellated visco foam contributes to an effective pressure ulcer prevention plan, offering immersion, comfort, and postural support




Internal air channels in the foam's surface help to reduce the buildup of heat and moisture on the skin

High frequency welded seams reduce the risk of fluid ingress, supporting infection control protocols

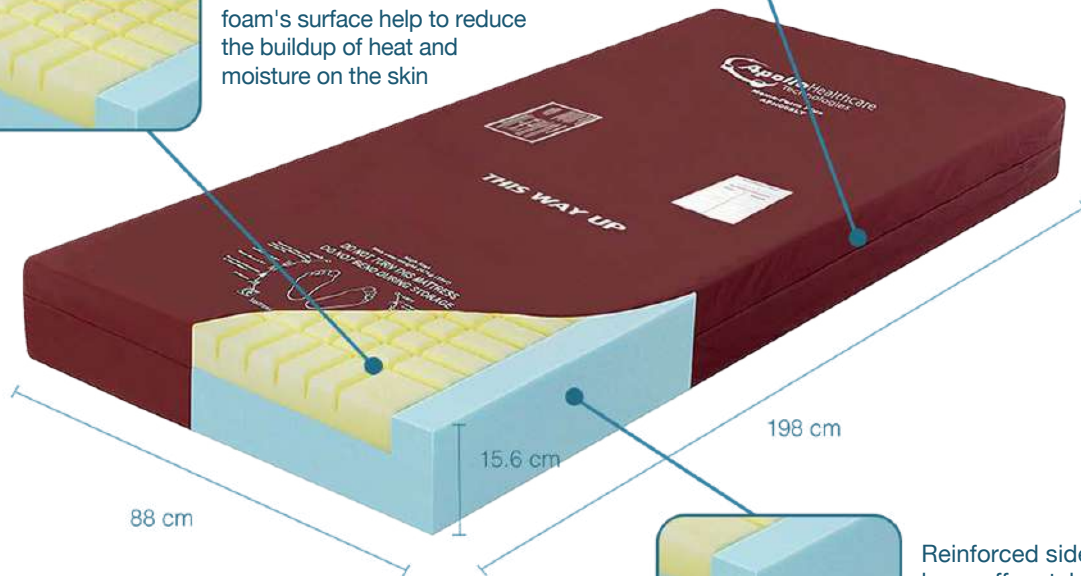

Reinforced side wall and base offer stability and support which assists with patient transfers

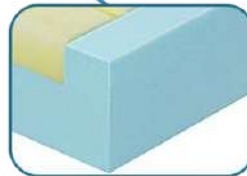

Specifications

Product Code - APH005LT

Dimensions: 198 x 88 x 15.6 cm

Mattress Weight: 16 kg

Maximum Weight Limit: 39 stone / 248 kg

Risk Factor: High Risk

Warranty: 2 years

Safety Standards: BS EN 597-1 Cig, BS EN 597-2 Match, BS 7177 Ignition Source 5

Due to ongoing research and development, Apollo Healthcare Technologies Ltd, reserve the right to change specifications without prior notice.

If you require more information, please call our sales team on Tel: +44 (0) 1924 614567



Apollo Healthcare Technologies Limited, Unit A1, Gildersome Spur, Leeds, West Yorkshire, LS27 7JZ
Tel: +44 (0) 1924 614567 Fax: +44 (0) 1924 607480 Email: sales@apollo-ht.co.uk www.apollo-ht.co.uk